


# NeoMagic<sup>®</sup>

## 1.9 FR Neonatal PICC

Neo-Magic<sup>®</sup> vascular access products feature a design focused on neonatal and pediatric patients. The redesigned 1.9 FR PICC gives clinicians working with neonatal patients a tool specialized to meet their needs.



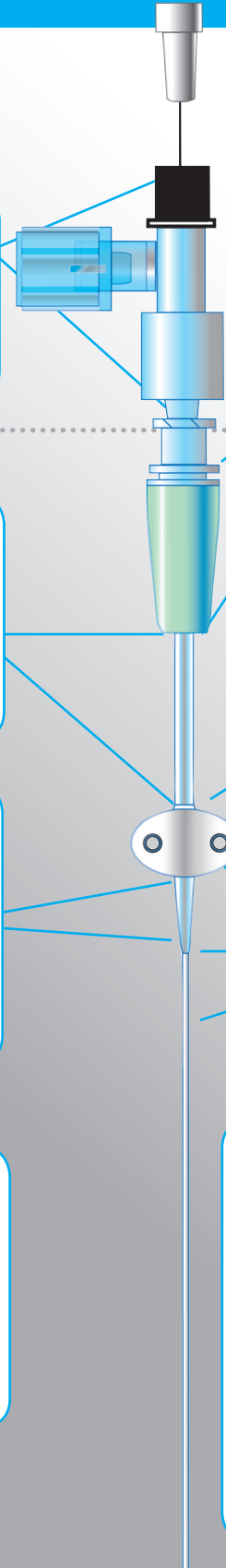

Patented Neo-Magic® EZ-Flush® stylet allows superior tip placement and flushing past vascular valves. This is the only stylet with this capability for very small catheters.

The Mini Profile Luer was designed with added security in mind. It is both sutured and sleeved and provides a non-bulky, comfortable profile.


The clear extension leg provides visibility of the catheter pathway even after securement. This feature allows for easier identification of potential problems such as blood retrograde and solution bubbling.

The flat-bottomed disc allows for the use of a taping technique made popular with the Neo PICC®. This technique is illustrated in the IFU. The holes in the disc to allow the use of the new StatLock® PICC II neonatal securement device.

The strain relief taper provides added strength in the most vulnerable area of the catheter. While this area still requires extra caution, the feature allows us to achieve the same specification strength results used for the catheter tubing.

The 3cm area from the strain relief to the marked skin exit point allows for use of popular securing techniques and provides a reference point for catheter tip movement.

The 1.9Fr silicone catheter tubing is numbered and marked with a very small yet visible pattern. This pattern limits surface imperfections that have been shown to cause fibrin attachment. The catheter is trimmable and available in lengths from 15 cm to 40 cm.

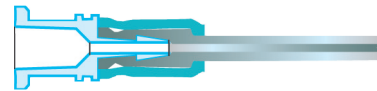

The unique internal hub lumen is tapered from the luer entrance to the exiting catheter lumen. This feature eases the insertion of a small diameter guide-wire in the event catheter repositioning is required.

Neo PICC® is a registered trademark of Arrow International. StatLock® is a registered trademark of C.R. Bard, Inc.

**Neo Medical Inc.**  
42514 Albrae St.  
Fremont, CA 94538  
(888)450-3334  
[www.neomedicalinc.com](http://www.neomedicalinc.com)